

Application Note



**Trace analysis of 29 perfluorinated compounds (PFAS)
in salmon trout using **AttractSPE® PFAS****

This application note describes an efficient solid phase extraction (SPE) method for the cleanup and trace analysis of 29 perfluorinated compounds (PFAS) in salmon trout. The method relies on **AttractSPE® PFAS**, a SPE cartridge, to catch perfluorinated compounds in complex matrices prior to their analysis by LC-MSMS.

Perfluorinated compounds (PFAS or PFCs) are a large family of molecules consisting of varying lengths of fluorocarbons chains with a functional group such as carboxylic or sulfonic acids attached. They have been widely used for more than 50 years in various products, such as firefighting foams, hydrophobic and nonstick coatings, or surfactants to cite a few examples. Their nature makes them particularly chemically inert and very resistant to degradation in environment. Some PFAS are classified as persistent organic pollutants (POPs) and are strongly associated with a variety of human disorders such as neurotoxicity, immune deficiency, and cancer[1]. In april 2023, the EU regulation commission 2023/915 set a limit of 0.2 – 2 µg/Kg for PFOA, PFOS, PFNA, PFHxS in fish muscle meat intended for the manufacture of baby food.

Solid phase extraction allows to purify and concentrate samples prior to analysis, which is particularly important to reach low concentrations or in complex matrices.

Table 1. List of the 29 tested PFAS

Compound	Abbreviation	CAS Number
Perfluorobutanoic acid	PFBA	375-22-4
Perfluoropentanoic acid	PFPeA	2706-90-3
Perfluorohexanoic acid	PFHxA	307-24-4
Perfluoroheptanoic acid	PFHpA	375-85-9
Perfluorooctanoic acid	PFOA	335-67-1
Perfluorononanoic acid	PFNA	375-95-1
Perfluorodecanoic acid	PFDA	335-76-2
Perfluoroundecanoic acid	PFUnA	2058-94-8
Perfluorododecanoic acid	PFDoA	307-55-1
Perfluorotridecanoic acid	PFTTrDA	72629-94-8
Perfluorotetradecanoic acid	PFTeDA	376-06-7
Perfluorohexadecanoic acid	PFHxDA	67905-19-5
Perfluorooctadecanoic acid	PFODA	16517-11-6
Perfluorobutanesulfonic acid	PFBS	375-73-5
Perfluoropentanesulfonic acid	PFPeS	2706-91-4
Perfluorohexanesulfonic acid	PFHxS	355-46-4

Perfluoroheptanesulfonic acid	PFHpS	375-92-8
Perfluorooctanesulfonic acid	PFOS	1763-23-1
Perfluorononanesulfonic acid	PFNS	68259-12-1
Perfluorodecanesulfonic acid	PFDS	335-77-3
Perfluoroundecanesulfonic acid	PFUdS	749786-16-1
Perfluorododecanesulfonic acid	PFDoS	79780-39-5
Perfluorotridecanesulfonic acid	PFTTrDS	791563-89-8
Hexafluoropropylene oxide dimer acid	HFPO-DA	13252-13-6
1H,1H, 2H, 2H-Perfluorooctane sulfonic acid	6:2 FTS	27619-97-2
N-ethyl perfluorooctanesulfonamidoacetic acid	N-EtFOSAA	2991-50-6
Perfluorooctanesulfonamide	FOSA	754-91-6
N-methyl perfluorooctanesulfonamide	MeFOSA	31506-32-8
4,8-Dioxa-3H-perfluorononanoic acid	ADONA	919005-14-4

1 Proceeding of the experiment

Sample preparation:

Four samples of salmon trout were spiked with 29 PFAS at 0.2-1µg/Kg (detailed concentrations table 3). A non spiked salmon trout sample was also extracted as a blank to measure the initial presence of PFAS.

1. Dice a salmon trout fillet, freeze overnight, then grind in a blender.
2. Weight 2g of ground salmon trout into a 15mL centrifuge tube. Add 10mL 0.05M KOH (in methanol), vortex and agitate slowly on a shaking table for 16h. Centrifuge at 2800rpm for 10 minutes. Collect the supernatant into another tube.
3. Add 10mL acetonitrile to the same trout sample, vortex and sonicate for 30 min. Centrifuge at 2800rpm for 10 minutes. Collect the supernatant and add it to the first supernatant fraction.
4. Add another 5mL 0.05M KOH (in methanol) to the same sample, vortex and shake by hand. Centrifuge at 2800rpm for 10 minutes. Collect the supernatant and add it to the first and second supernatant fraction.
5. Add 1mL ultrapure water to the obtained solution.
6. Concentrate the solution obtained to around 2.5mL using nitrogen evaporation at 55°C. Adjust to ~48mL with ultrapure water and vortex. Adjust pH at 6-7 with formic acid or 3% NH₃ to form the loading solution prior to SPE.

**PURIFICATION WITH A 6mL
ATTRACTSPE® PFAS - 200mg CARTRIDGE WITH
60mL RESERVOIR AND ADAPTER**

Conditioning/Equilibration

1. 15mL 1% NH₃ in methanol
2. 5mL 0.3M formic acid in water

Loading

1. ~48mL of loading solution

Washing

1. *2x5mL ultrapure water
2. *5mL 0.05M formic acid in 50/50 water/methanol (v/v)
3. Dry cartridge for 15 seconds by applying full vacuum

Elution

1. *5mL 1% NH₃ in methanol

Analysis

Add 25µL of acetic acid to the elution and homogenize prior LC-MS/MS analysis

***Prior to washing/elution, rinse the walls of sample container with the solution then pass through the cartridge. For the elution step, also rinse the walls of the reservoir prior passing through the cartridge.**



Some molecules of interest are ubiquitous and so likely to be found in mobile phase or LC parts, leading to contamination during LC-MS/MS analysis. Fitting **SilactHPLC Delay – PFAS** as a delay column upstream of the injector (Figure 1) shifts the contaminant retention time, allowing the analyte of interest to be measured in the sample without contamination (Figure 2). This solution, in comparison with the monitoring of contamination from the LC device, is very easy to implement and is cost effective.

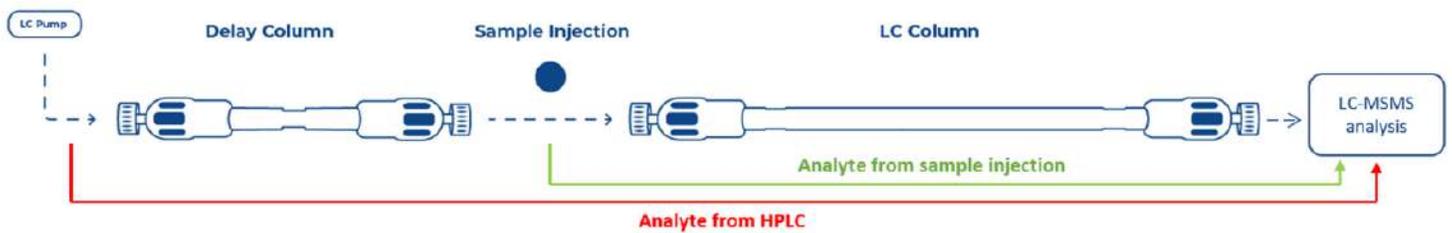


Figure 1. Installation of Delay column for LC analysis

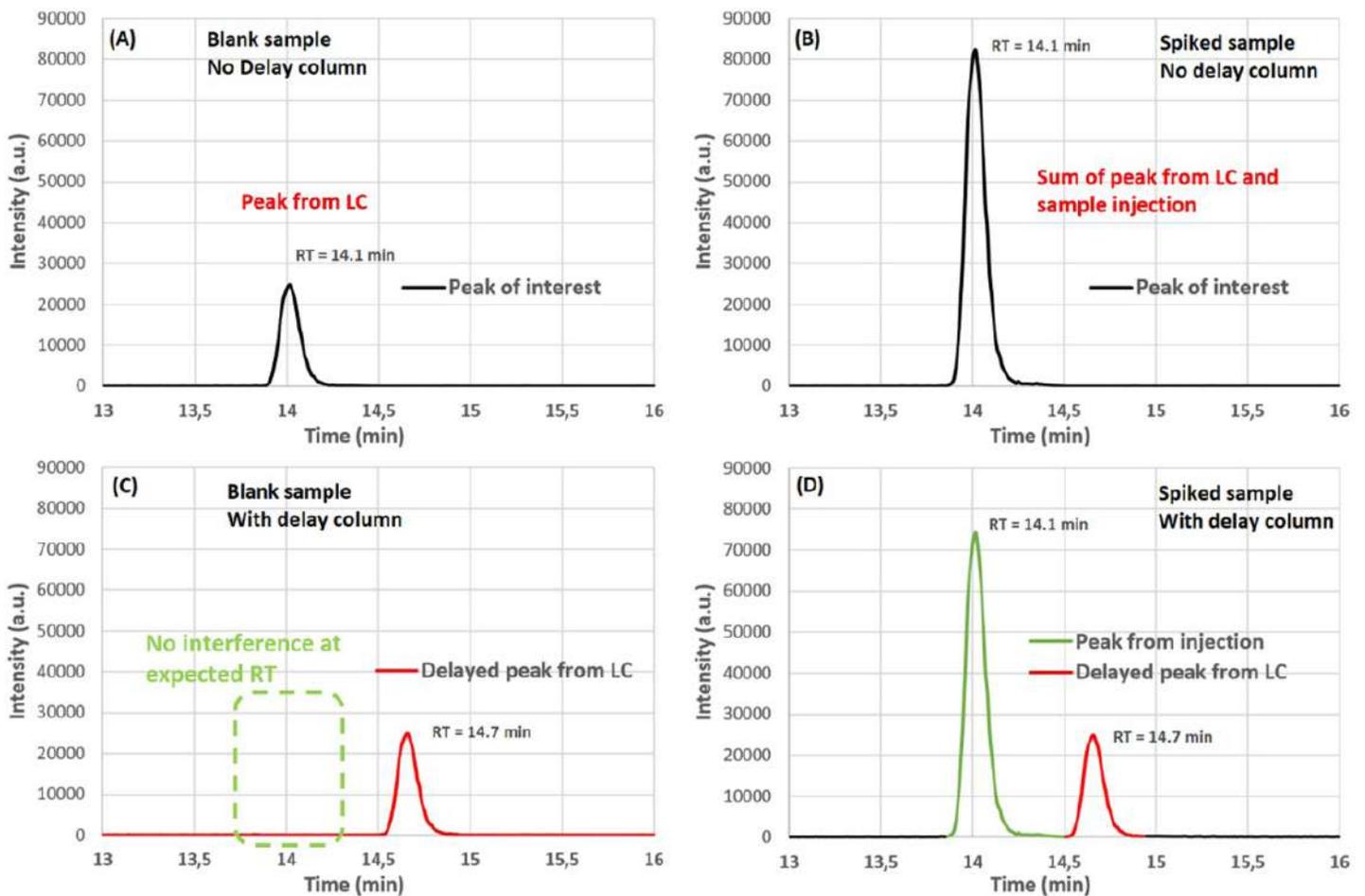


Figure 2. Effects of a delay column on samples. A: Blank sample without delay column. B: Spiked sample without delay column. C: Blank sample with delay column. D: Spiked sample with delay column.

Table 2. LC-MS/MS conditions of analysis

LC Conditions			MS/MS Conditions				
LC Dionex U3000			Qtrap 4000 ESI- MS/MS				
Column : Siliachrom dtC18 150*2.1mm at 40°C Delay column : SilactHPLC Delay-PFAS 50*2.1mm			Curtain gas : 30				
Injection volume : 5µL			CAD : High				
T° sampler : 10°C			IS : -4500V				
Flow rate : 0.25mL/min			Temperature : 550°C				
			GS1/GS2 : 50/50				
Time (min)	10mM Ammonium Acetate (in water)	Methanol	Analyte	Retention time (min)	Q1 (m/z)	Q3 (m/z)	CE (V)
0	60%	40%	PFBA	3.2	213.0	168.8	-14
			PFPeA	4.9	263.0	218.8	-12
1	60%	40%	PFHxA	10.1	313.0	268.9	-14
			PFHpA	13.0	363.0	318.8	-16
20	10%	90%	PFOA	15.0	413.1	368.9	-14
			PFNA	16.7	463.0	418.9	-16
30	10%	90%	PFDA	18.1	513.0	469.0	-18
			PFUnA	19.3	563.1	519.0	-16
31	60%	40%	PFDoA	20.3	613.1	569.1	-18
			PFTTrDA	21.1	663.1	619.0	-18
35	60%	40%	PFTeA	22.0	712.9	668.9	-18
			PFHxDA	23.2	813.0	769.0	-20
			PFODA	24.3	913.0	869.0	-22
			PFBS	6.2	299.0	79.8	-52
			PFPeS	10.6	349.0	79.9	-68
			PFHxS	13.2	399.0	79.9	-74
			PFHpS	15.1	448.9	80.0	-90
			PFOS	16.7	499.0	80.1	-84
			PFNS	18.1	549.0	80.0	-94
			PFDS	19.2	599.0	79.8	-110
			PFUdS	20.3	649.0	79.9	-116
			PFDoS	21.1	699.0	79.8	-130
			PFTTrDS	21.9	749.0	80.0	-130
			HFPO-DA	11.1	285.1	168.7	-12
			6:2 FTS	14.9	427.1	406.8	-34
			N-EtFOSAA	19.4	584.1	418.8	-30
			FOSA	19.7	498.1	77.8	-70
			MeFOSA	22.1	512.1	168.8	-38
			ADONA	13.3	377.0	250.9	-18

3 Results

The analytes were simultaneously analysed by LC-MS/MS. The results obtained are presented in the table below.

Table 3. Recovery obtained for tested analytes (ND : not detected)

Compound	Salmon trout spiked at (µg/kg)	Salmon trout		
		[C] in blank (µg/kg)	% Recovery	% RSDr (n=4)
PFBA	0.20	ND	95%	6%
PFPeA	0.20	ND	116%	5%
PFHxA	0.20	ND	103%	3%
PFHpA	0.20	ND	102%	4%
PFOA	0.20	ND	102%	5%
PFNA	0.20	ND	97%	13%
PFDA	0.20	ND	98%	6%
PFUnA	0.20	ND	87%	4%
PFDoA	0.20	ND	87%	3%
PFTTrDA	0.20	0.08	96%	17%
PFTeA	0.20	ND	96%	11%
PFHxDA	0.20	ND	87%	10%
PFODA	0.20	ND	88%	12%
PFBS	0.20	ND	94%	7%
PFPeS	0.20	ND	93%	4%
PFHxS	0.20	ND	97%	6%
PFHpS	0.20	ND	92%	4%
PFOS	0.20	0.17	89%	9%
PFNS	0.20	ND	101%	8%
PFDS	0.20	ND	94%	5%
PFUdS	0.20	ND	90%	6%
PFDoS	0.20	ND	93%	6%
PFTTrDS	0.20	ND	95%	4%
HFPO-DA	0.20	ND	96%	10%
6:2 FTS	1.00	0.02	88%	11%
N-EtFOSAA	0.40	ND	105%	13%
FOSA	0.20	0.05	77%	10%
MeFOSA	0.40	ND	84%	10%
ADONA	0.20	ND	98%	6%

CONCLUSION

AttractSPE® PFAS has been successfully used for the enrichment and cleanup of 29 perfluorinated compounds in salmon trout. The method has shown excellent performances with recoveries from **77% to 116%** and a very good repeatability.

Furthermore, the use of **SilactHPLC DELAY - PFAS** as delay column allows to avoid any PFAS interference during LC-MS/MS analysis. Particular attention must be paid to check that the laboratory environment does not contaminate samples and lead to false positives. Some simple precautionary steps are described in the application note (e.g., the use of a delay column). For routine analysis, the use of internal standards to correct the potential matrix effects and adsorption of the largest PFAS is advised.

References

1. **Impact of Perfluorinated Compounds on Human Health, 2014 Academy for Environment and Life Sciences.**

Product references:




***AttractSPE® PFAS**
6mL - 200mg
Quantity : 50/pk
PFAS-50.S.6.200



***AttractSPE® PFAS GCB-WAX**
6mL - 50mg/200mg
Quantity : 50/pk
PFAS-GW-50.S.6.50.200



SilactHPLC DELAY - PFAS
50*2.1mm (5µm)
for PFAS analysis
Quantity : 1 unit
DELAY-PFAS-50.2.1

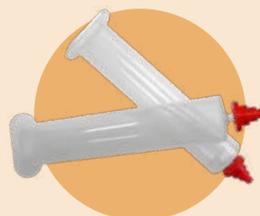

**kit of 12 reservoirs 60ml and
adapters for use with 1,3 & 6 mL
columns**
Quantity : 1 unit
ACC-ARI



SPE Vacuum Manifold
12-port model
Quantity : 1 unit
ACC-MAN2

***Other formats and sorbent volumes available
in our catalogue**